**PHYSIO RADIANCE**

P a r i s

*VISAGE+*

The evolution of beauty

USER MANUAL

# PHYSIO RADIANCE

P a r i s

Congratulations on your purchase of the Physio RadianceVisage+, your ultimate skin care solution. Physio RadianceVisage+ will help you look your best with beautiful and flawless skin. Let Physio RadianceVisage+ be your new transformation. Start seeing results today.

---

## VISAGE+ TABLE OF CONTENTS

1. SAFETY GUIDELINES	<b>3</b>
2. ASSEMBLING	<b>4</b>
3. GET STARTED	<b>5</b>
• TECHNOLOGY USED	
4. SPECIFICATIONS	<b>5</b>
5. USING THE PHYSIO RADIANCEVISAGE+	<b>6</b>
• GETTING TO KNOW THE DEVICE	
• USING CONTROLLER	
• CHARGING THE DEVICE	
• USING THE DEVICE	
• AFTER USAGE CARE	
6. FAQs	<b>11</b>
7. PRODUCT WARRANTY	<b>12</b>

# SAFETY GUIDELINES

- Always follow the basic precautions when using this device.
- Do not use the device if you have any electronic implants such as a cardiac pacemaker or defibrillator as the device may interfere with the proper functioning of these electronic implants.
- Do not use the device if you are connected to any other electrical equipment with high frequency or diathermy to avoid the risk of any malfunctions due to electrical interference.
- Do not use the device over areas with braces, fillings or metal implants as you may experience some discomfort or a tingling sensation.
- Do not use or place the device directly over your eyes.
- Do not use the device over swollen, infected or inflamed areas and skin eruptions.
- Do not use over broken skin, lacerations, skin disorders and infectious diseases.
- Do not use the device if you suffer from serious illness or terminal health conditions.
- Do not use the device if you have numbness or insensitivity to electrical or heat treatments.
- Do not use the device if you are under medical supervision for any cognitive dysfunction.
- Do not use on children.
- Keep out of reach of children.
- Do not use on hyper apprehensive or highly strung users.
- Ensure that there is ample amount of skin care applied to the face and neck prior to using the device.
- Do not use if you are pregnant.
- Do not use the device while charging.
- Do not use when hands are wet.
- Do not submerge the device in water.
- Consult your physician before using the device.

# ASSEMBLING

The following items are provided in your Physio RadianceVisage+ kit:

1. Physio RadianceVisage+ Device
2. Physio RadianceVisage+ Charging Cradle
3. Physio RadianceVisage+ Travel Pouch
4. Universal Power Adaptor & USB Cable
5. User manual

# GET STARTED

## TECHNOLOGY USED

Physio RadianceVisage+ is a unique facial device equipped with a combination of Galvanic technology, Radio Frequency and Chromo Therapy.

This simple and easy-to-use facial device is designed for an array of facial and beauty treatments, helping to cleanse the skin, introduce active substances to the skin and enhance regeneration of elastin fibers and collagen tissues. Physio RadianceVisage+ has a user friendly control panel and an LED display panel with preset functions and intensity levels to suit all skin types. The device is mounted with four Titanium Facial Pads and three coloured LED Lights.

**Galvanic technology** (C, M) offers an ultimate facial care procedure through Desincrustation and iontophoresis. Desincrustation technology deeply cleanses and moisturises the skin while Iontophoresis technology rejuvenates the skin by infusing active substance into the deeper layers of skin tissue.

**Radio Frequency** (R, L) offers facial contour lifting by re-programming the regeneration of the collagen tissues and muscle fibres, hence reducing and smoothing fine lines and wrinkles on the face.


**Chromo Therapy** (LED) is a recently developed facial treatment using three LED lights – red, yellow and blue – set at specific wavelengths to suit different types of skin levels.

## SPECIFICATIONS


<b>Product name</b>	: Physio RadianceVisage+
<b>Input Voltage</b>	: 100-250 V
<b>Output Voltage</b>	: 5 V
<b>Battery</b>	: Rechargeable Lithium Battery
<b>Output Current</b>	: 1,000 mA
<b>Dimensions</b>	
<b>Length</b>	: 18 cm
<b>Width</b>	: 4.5 cm
<b>Height</b>	: 3.3 cm
<b>Weight</b>	: 118 g

# USING THE PHYSIO RADIANCE VISAGE+


## GETTING TO KNOW THE DEVICE

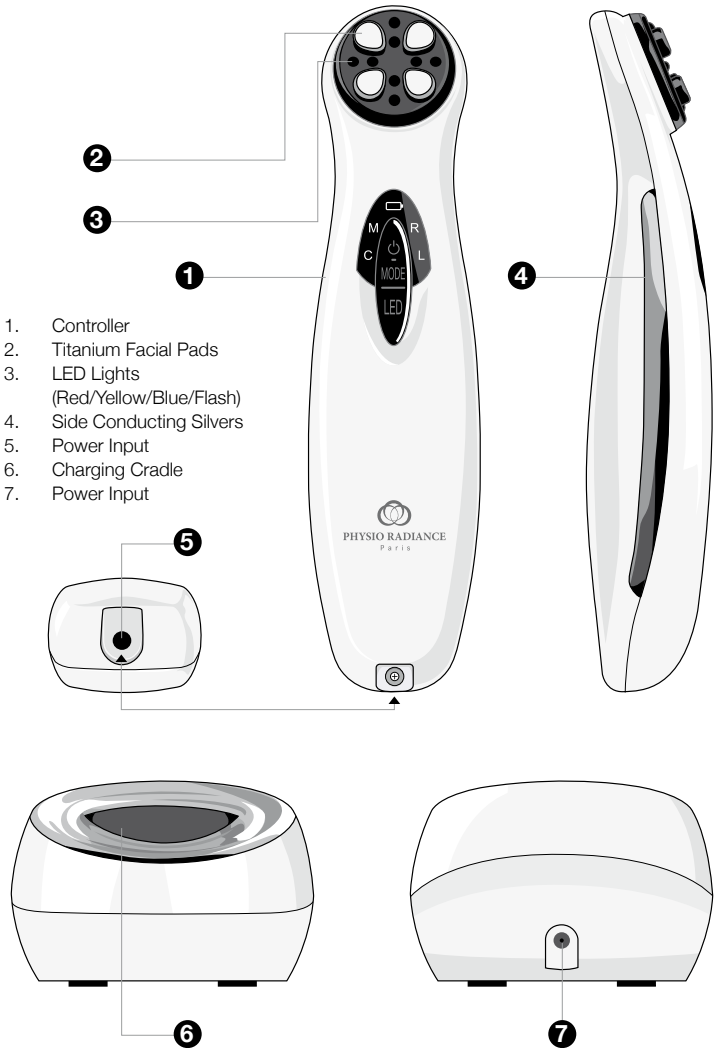
## USING CONTROLLER

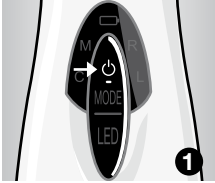
**On/Off**

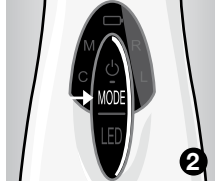
**Mode selection**  
(Galvanic/Radio Frequency)

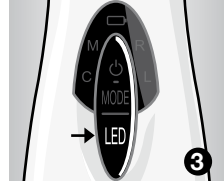
**LED Lights selection**

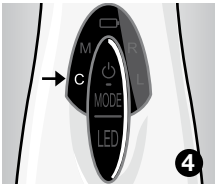
**C Cleansing**

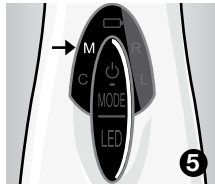
**M Moisturising**

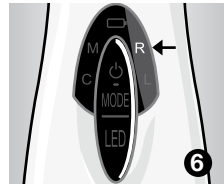**R Radiance**

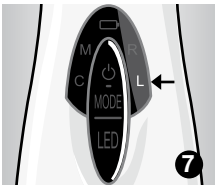**L Lifting**

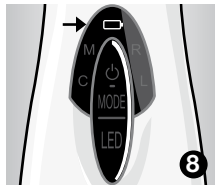**Battery Indicator**

## CHARGING THE DEVICE

1. Ensure all accessories and components are complete.
2. Charge the device with the provided adaptor only.
3. To charge, insert the device into the Cradle and make sure the device sits snugly in the cradle. Once the adaptor is plugged into the socket, insert the jack into the power input situated at the back of the cradle. While charging, the battery indicator light will turn red indicating that the charging is in progress.
4. Once fully charged, the indicator light will turn green.

**Note:** It will normally take 1.5-2 hours to charge a fully drained battery.


**DO NOT USE THE DEVICE WHILE CHARGING.**

# USING THE DEVICE

## BEFORE USING THE DEVICE


- Prior to using Physio RadianceVisage+, ensure that the skin on and around the face and neck area is clean for optimal result.
  - The duration of each Mode is three (3) minutes. The usage duration time will be indicated by the beeping sounds, first beep being the one (1) minute mark, the second beep being the two (2) minute mark, and the last beep being the three (3) minute mark. The device will automatically turn off after the final beep
  - During every procedure, ensure that all four titanium facial pads are in full contact with the skin at all times.
  - Ensure you touch the side conducting silvers on both side of the device at all times to complete the circuit.
  - After every session, wipe or clean excess gel/skin care products off the device and titanium pads with a clean damp cloth to remove traces of residues.
- 

## STEP 1 DEEP CLEANSING (C)

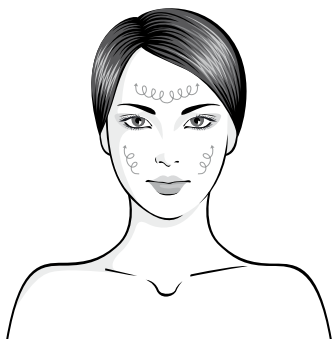- Spread an ample amount of Physio Radiance Visage+ gel onto the face and neck area.
- Turn on Physio RadianceVisage+ and set to C (Cleansing) with Red LED light. Gently glide the device around your face and neck in small rotating motions (You can choose to use Blue LED light if you have congested skin or acne).
- Start from the chin then move to one side of the cheek then up to the nose area, followed by the forehead area and down to the other side of the cheek. Continue with the procedure until you are back to the chin. Repeat the procedure for three (3) minutes.
- It is mandatory to rinse your face with warm water and pat dry with facial tissue or cloth before proceeding to the next step.

## STEP 2 MOISTURISING (M)

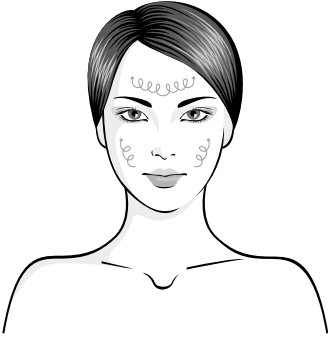
- Spread an ample amount of serum/moisturiser onto the face and neck area.
- Turn on Physio RadianceVisage+ and set to M (Moisturising) with Red LED light. Gently massage your face and neck in rotating motions.
- Start from the chin then move to one side of the cheek then up to the nose area then up to the forehead area and down to the other side of the cheek. Continue with the procedure until you are back to the chin. Repeat the procedure for three (3) minutes.

## STEP 3 RADIANCE (R)

- Spread an ample amount of serum/moisturiser onto the face and neck area.
- Turn on Physio RadianceVisage+ and set to R (Radiance) with Red LED light. Gently massage your face and neck in rotating circles.
- Start from the chin then move to one side of the cheek then up to the nose area then up to the forehead area and down to the other side of the cheek. Continue with the procedure until you are back to the chin. Repeat the procedure for three (3) minutes.

## STEP 4 LIFTING (L)

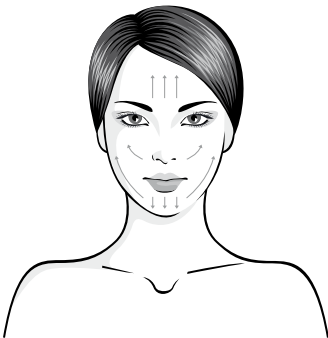- Spread an ample amount of serum/moisturiser onto the face and neck area.
- Turn on Physio RadianceVisage+ and set to L (Lifting) with Red LED light.
- Start from the lower lip area down towards the chin then begin at the corner of the mouth from jawline up towards the temple area. Gently massage in an upward motion. Repeat the procedure on the other side of the face for three (3) minutes.

## RECOMMENDED PROGRAMS

<b>SKIN TYPE/ CONDITION</b>	<b>CLEANSING (C)</b>	<b>MOISTURISING (M)</b>	<b>RADIANCE (R)</b>	<b>LIFTING (L)</b>	<b>LIGHT</b>
Normal	X	X	X	X	Red
Dry		X	X	X	Red
Oily & Congested	X	X	X	X	Red
Matured	X	X	X	X	Red
Dehydrated	X	X	X	X	Red
Acne	X				Blue
Sensitive		X	X	X	Yellow
Dark Eye Rings	X	X	X	X	Flash

## AFTER USAGE CARE

- Wipe or clean excess gel/skin care products off the device and the titanium pads with a clean damp cloth to remove any traces of residue after every session.
- Store the device and its components in a dry and dust-free environment.
- Store in temperatures between -10 °C to 35 °C.
- Do not spill any liquid on the device, controller, power adaptor, cradle or any electrical wires.
- Do not use in wet conditions e.g. rain or in a shower.
- Do not clean the device with detergents or chemicals of any sort.
- Do not submerge the device in water.

# FAQs

## **How does a Radio Frequency Facial work?**

Radio Frequency technology works through radio wave currents that penetrate deep into the dermis (deepest layer of the skin) and internal layers of the skin, stimulating the tightening of the underlying tissues from which fine lines, wrinkles, and bad skin texture originate as part of the natural ageing process.

## **Does the procedure hurt?**

User experience may vary due to different levels of pain tolerance in individuals. The radio wave currents together with the LED light therapy projects heat onto the skin and may cause minimal discomfort.

## **Some parts of my face and neck turned red after use, is this normal?**

Due to radio wave currents and heat, users may experience mild redness. This should go away in 1-2 hours. If redness persists, please consult your physician immediately.

## **How long does the effects/results last?**

Results vary from user to user. They can last up to 24 months or more based on the user's natural ageing process.

## **How many times a week should I use this device?**

For best results, use 7 – 14 times a week.

## **How long is each session?**

A full session will take twelve (12) minutes while each individual session takes three (3) minutes each.

## **Who can use this device?**

This device is best for those in their mid-20s and 30s-50. Users in their 60s and above may require more frequent procedures.

## **What kind of skin type is best suited for the device?**

The device can be used with all skin types.

## **What kind of skin care should I use with the device?**

The device can be used with any skin care products that contain water-based ingredients. For best and optimal results, we recommend you to use the device with Physio Radiance skin care products.

# PRODUCT WARRANTY

To be covered by this Limited Warranty, the product must have been purchased directly from QNET Limited ('QNET'), or from a registered Independent Representative (IR). This Limited Warranty extends only to the original purchaser and is not assignable or transferable.

QNET WARRANTS TO YOU, THAT UNDER NORMAL USE, THE PRODUCT WILL BE FREE FROM DEFECTS IN MATERIAL AND WORKMANSHIP FOR A PERIOD OF ONE (1) YEAR FROM THE DATE OF PURCHASE, ('WARRANTY PERIOD'). IF YOU DISCOVER A DEFECT COVERED BY THIS LIMITED WARRANTY DURING THE WARRANTY PERIOD, AND RETURN THE PRODUCT TO US, WE WILL REPAIR OR REPLACE THE PRODUCT, AT OUR SOLE OPTION, USING EITHER NEW OR REMANUFACTURED COMPONENTS. THERE IS NO WARRANTY OR CONDITION OF ANY KIND AS TO ANY DEFECTS DISCOVERED AFTER THE WARRANTY PERIOD EXPIRES.

This Limited Warranty covers defects in manufacturing that arise from the correct use of the product per the product's manual and instructions. It does not cover damage caused by abuse, misuse, accidents, modification or repair by any persons other than QNET Ltd.

THIS LIMITED WARRANTY IS THE ONLY WARRANTY FOR THE PRODUCT AND THERE ARE NO OTHER EXPRESS WARRANTIES, ORAL OR WRITTEN, PROVIDED BY QNET. THIS LIMITED WARRANTY DOES NOT GUARANTEE SPECIFIC SKIN IMPROVEMENT.

## **LIMITATION OF REMEDIES**

YOUR EXCLUSIVE REMEDY, IN LIEU OF ALL INDIRECT, INCIDENTAL, CONSEQUENTIAL, EXEMPLARY, SPECIAL OR OTHER DAMAGES OF ANY KIND, INCLUDING NEGLIGENCE, IS LIMITED TO THE REPAIR OR REPLACEMENT OF THE PRODUCT, IF DEFECTIVE, UNDER THE TERMS AND CONDITIONS OF THIS LIMITED WARRANTY.

QNET shall bear no other expenses.

QNET reserves the right to refuse to repair or replace the product if it is returned in an unsanitary condition, and reserves the right to substitute a product of equal or better quality if an identical product is not available at the time of replacement.

## **GOVERNING LAW**

The laws of Hong Kong shall apply to this Limited Warranty.

## **HOW TO GET SERVICE UNDER THIS LIMITED WARRANTY**

To obtain service under this Limited Warranty, please contact the Global Support Centre and keep the product serial number for tracking purposes.

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Copyright © QNet Ltd / 2014. All rights reserved.
Copyrights and/or other intellectual property rights on all designs, graphics, logos, images, photos, texts, trade names, trademarks, etc., in this publication are reserved. The reproduction, transmission, or modification of any part of the contents of this publication is strictly prohibited.
**QNet Ltd**